

OLYMPUS

Your Vision, Our Future

26 INCH HIGH-DEFINITION LCD MONITOR

OEV262H


The Ideal Synergy with Your Olympus Endoscopy System.



OE261H

Main Features

- 26-inch full-HD LCD panel with a 16:9 aspect ratio, high brightness, high contrast, and high gradation images.
- FLIP function (such as mirror and 180° rotation) which provides a suitable display pattern and monitor layout for the procedure.
- Multimodality display capability, including picture-in-picture (PIP), picture-out-picture (POP), and clone out*¹ with various image-size combinations.
- Various input/output terminals, including 3G/HD/SD SDI (2x), DVI-I (2x), HD15, Y/C, and Composite.
- Eco-friendly specifications thanks to low power consumption, various power-saving modes, and a lightweight and thin body compared to its predecessor*².
- A.I.M.E (Advanced Image Multiple Enhancer)*³ which provides support by enhancing the image quality of the OLYMPUS endoscope system.


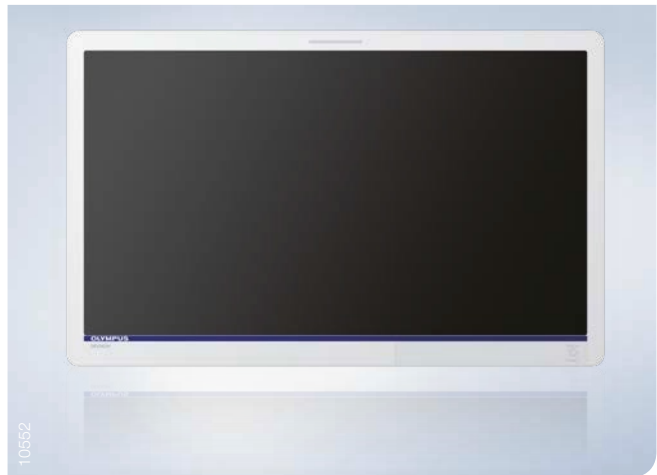


- * 1 Clone out is a function outputting a HD video signal (1080i/1080p) as displayed on the monitor, including the PIP/POP image to the second monitor and recording device.
- * 2 HIGH-DEFINITION LCD MONITOR OE261H.
- * 3 A.I.M.E. (Advanced Image Multiple Enhancer) is a function that significantly improves the resolution of fine structures without adding additional noise to the image.

Specifications

LCD Panel	Screen size	26 inches
	Display devices	a-Si TFT active matrix
	Resolution	1,920 × 1,080 pixels (full HD)
	Contrast ratio	1,400:1
	Aspect ratio	16:9
	Viewing angle	178° (horizontal and vertical)
	Number of colors	1.07 billion
Size	Dimensions	626.8 (W) × 395.2 (H) × 79.6 (D) mm
	Weight	7.8 kg
Power Supply	LCD monitor	AC IN: 100–240 V, 50/60 Hz, 0.96–0.47 A DC IN: 24 V, 3.6 A, 5 V, 0.03 A (supplied by optional AC adapter)
	Input	HD/SD SDI (1)
SDI (2)		BNC connector
VIDEO		BNC connector
Y/C		4-pin mini-DIN
HD15		D-sub 15-pin
DVI (1)		DVI-I connector
DVI (2)		DVI-I connector
Aux in (HD/SD-SDI)		BNC connector
Output	HD/SD SDI (1)	BNC connector with active-through output
	SDI (1) SDI (2)	BNC connector with active-through output
	VIDEO	BNC connector with loop-through output and auto 75 termination
	Y/C	4-pin mini-DIN with active-through output
	HD15	D-sub 15-pin with active-through output
	DVI (1)	DVI-I connector with active-through output
	DVI (2)	DVI-I connector with active-through output
	Clone Out (3G/HD-SDI)	BNC connector
Remote	RS-232C connector	D-sub 9-pin
	GPI connector	RJ-45

Dimensions (mm)

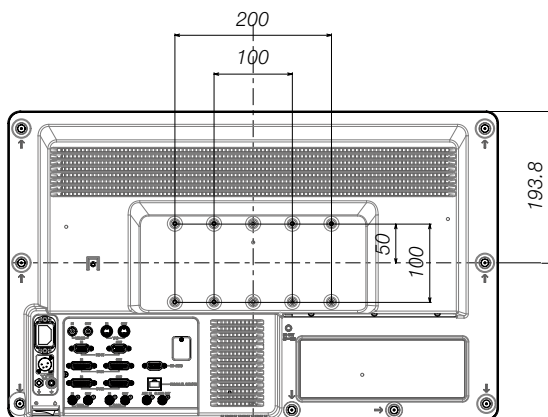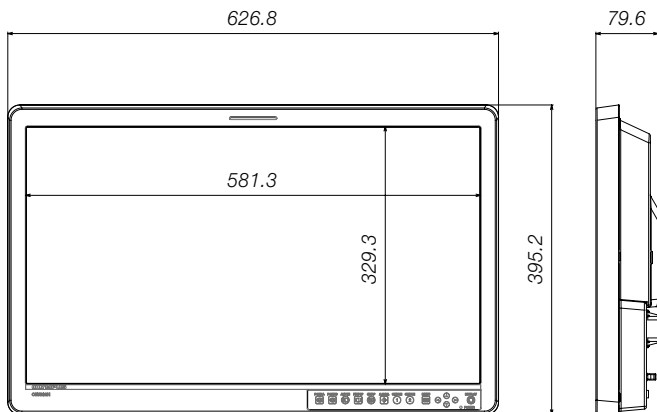

Specifications, design, and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS

OLYMPUS EUROPA SE & CO. KG
 Postbox 10 49 08, 20034 Hamburg, Germany
 Wendenstrasse 14–18, 20097 Hamburg, Germany
 Phone: +49 40 23773-0, Fax: +49 40 233765
 www.olympus-europa.com